



PLANETARY BALL MILL PM 200

The Planetary Ball Mill PM 200 is a powerful benchtop model with 2 grinding stations for grinding jars with a nominal volume of 12 ml to 125 ml.

The extremely high centrifugal forces of Planetary Ball Mills result in very high pulverization energy and therefore short grinding times.


The PM 200 can be found in virtually all industries where the quality control process places the highest demands on purity, speed, fineness and reproducibility.

The mill is ideally suited for tasks in research like mechanochemistry (mechano-synthesis, mechanical alloying and mechanocatalysis), or ultrafine colloidal grinding on a nanometer scale, as well as for routine tasks such as mixing and homogenizing soft, hard, brittle or fibrous materials.



TWO GRINDING STATIONS FOR STANDARD APPLICATIONS

- | Max. speed 650 rpm
- | Up to 10 mm feed size and 0.1 µm final fineness
- | 2 grinding stations for jars from 12 ml up to 125 ml, jars of 12 and 25 ml can be stacked (two jars each)
- | GrindControl to measure temperature and pressure inside the jar.
- | Aeration lids to control the atmosphere inside the jar
- | Storable SOPs and cycle programs, 5 different jar materials for dry and wet grinding

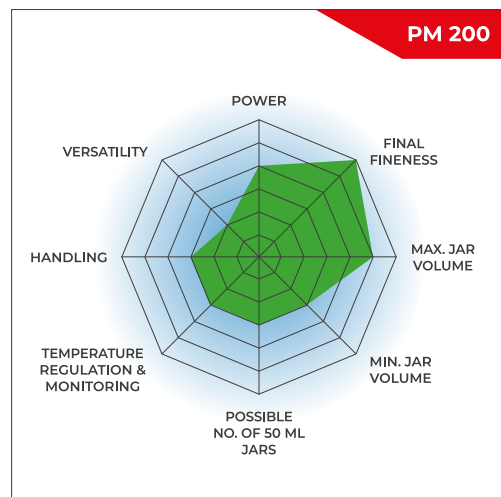

FAST & POWERFUL


- | Loss-free size reduction down to the submicron range
- | Wet grinding yields particle sizes in the nanometer range (<100 nm)
- | Variable speed from 100 to 650 rpm, speed ratio 1:-2
- | Grinding with up to 33.3 x acceleration of gravity
- | Batch-wise processing with max. 2 x 50 ml sample
- | Wide range of materials for contamination free grinding

PM 200

PLANETARY BALL MILL PM 200

REPRODUCIBILITY, SAFETY AND EASY HANDLING

- | Reproducible results due to speed control
- | Easy and safe clamping of grinding jars
- | The Safety Slider prevents starting the machine without securely clamped jars
- | Perfect stability on the lab bench thanks to FFCS technology
- | Innovative counterweight and imbalance sensor for unsupervised operation
- | Comfortable parameter setting via display and ergonomic 1-button operation
- | Automatic grinding chamber ventilation
- | 10 SOPs can be stored, programmable starting time
- | Power failure backup ensures storage of remaining processing time

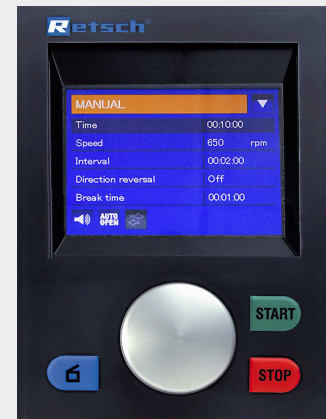

SETTINGS & OPTIONS

- | Dry and wet grinding possible
- | Suitable for long-term trials, 99 h max.
- | Interval operation allows for cooling breaks
- | Direction reversal helps to minimize caking effects

THE BEST
ALTERNATIVE TO A
RETSCH PLANETARY
BALL MILL? A
RETSCH MIXER MILL.



Benefit from particularly ergonomic handling while achieving the same finenesses down to the nanometer range.

PLANETARY BALL MILL PM 200

SAFETY FIRST: JAR CLAMPING

Operation of the RETSCH planetary ball mills is particularly safe. They feature a robust Safety Slider which ensures that the mill can only be started after the grinding jar has been securely fixed with a clamping device. The self-acting lock ensures that the jar is seated correctly and securely. This proven solid mechanical system is less failure-prone than electronic solutions - the user has full access to the sample at any time. When the electronic system fails, it is not possible to unlock the jars, for example.



[Click to view video](#)

PLANETARY BALL MILL PM 200

WET AND NANO-SCALE GRINDING WITH THE PM 200

Wet grinding is used to obtain particle sizes below 5 μm , as small particles tend to get charged on their surfaces and agglomerate, which makes further grinding in dry mode difficult. By adding a liquid or dispersant the particles can be kept separated.

To produce very fine particles of 100 nm or less (nano-scale grinding) by wet grinding, friction rather than impact is required. This is achieved by using a large number of small grinding balls which have a large surface and many friction points. The ideal filling level of the jar should consist of 60 % small grinding balls.


For more details on jar filling, wet grinding and sample recovery watch the video.



[Click to view video](#)

The video shows wet grinding in the Planetary Ball Mill PM 100.

The graphic shows the result of grinding barium titanate at 500 rpm in the PM 200. After 5 h of pulverization in heptane and oleic acid mixture with 0.5 mm grinding balls, the D90 value of the original sample was reduced from 15 μm to 95 nm.

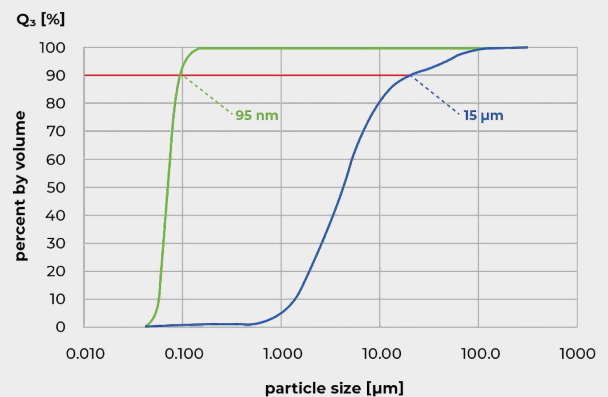

Grinding of barium titanate in heptane and oleic acid mixture with 0.5 mm grinding balls.

Blue curve: original sample; green curve: pulverized sample after 5 h.

PLANETARY BALL MILL PM 200

EASYFIT GRINDING JARS FOR EXCELLENT RESULTS

The performance and the result of sample preparation are also determined by the choice of the grinding jar and its ball charge. The EasyFit range of jars has been specially designed for extreme working conditions such as long-term trials, even at maximum speed of 800 rpm, wet grinding, high mechanical loads and maximum speeds as well as for mechanical alloying. This line of jars is suitable for all RETSCH planetary ball mills.

The new EasyFit grinding jar series features a structure on the bottom of the 50-500 ml jars called Advanced Anti-Twist (AAT). This ensures that the jars are tightly fixed without the risk of twisting, even at high speed, and that wear and tear is drastically reduced. Secure clamping of the jars is made much easier: to find the correct clamping position, a maximum twist of 60° is required.

The geometry of the EasyFit jars in the 50 ml and 250 ml sizes has been enlarged in diameter and reduced in height compared to the previous "comfort" models. This offers two advantages: better grinding results and interchangeable lids, as there are only three diameter dimensions for the entire grinding jar range.

Diameter categories

- | Diameter 1: 12 ml and 25 ml grinding jars
- | Diameter 2: 50 ml, 80 ml and 125 ml grinding jars

- | Available jar sizes: 12 ml / 25 ml / 50 ml / 80 ml / 125 ml
- | Innovative Advanced Anti-Twist (AAT) function ensures secure fit of grinding jars
- | High flexibility thanks to suitability of three lid sizes for all seven jar sizes
- | Pressure-tight and dust-proof O-ring sealing prevents material spillage
- | Jars and balls available in 5 materials: hardened stainless steel, tungsten carbide, agate, sintered aluminium oxide, zirconium oxide
- | Stainless steel protective jacket for agate, sintered aluminum oxide, zirconium oxide and tungsten carbide grinding jars
- | A groove between jar body and lid allows for easy opening of the lid, e. g. with the help of a spatula, if there are underpressure effects inside the jar



JARS & LIDS FOR SPECIAL APPLICATIONS

- | For colloidal or wet grinding, the use of a grinding jar with a special closure device is recommended
- | The special closure device is designed for ergonomic handling
- | Aeration lids are designed for working under inert atmosphere, for example if oxygen can influence the grinding process or the mechanosynthesis. The lids allow the introduction of gases like argon or nitrogen into the grinding jar.



Aeration lid



[Click to view video](#)

Video: Aeration lid

PLANETARY BALL MILL PM 200

RECOMMENDED JAR FILLINGS

To produce optimum grinding results, the jar size should be adapted to the sample amount to be processed. The grinding balls are ideally sized 3 times bigger than the largest sample piece. Following this rule of thumb, the number of grinding balls for each ball size and jar volume is indicated in the table below. To pulverize, for example, 50 ml of a sample consisting of 3 mm particles, a 125 ml jar and grinding balls sized at least 10 mm or larger are recommended. According to the table, 30 grinding balls are required.

Grinding jar nominal volume	Sample amount	Max. feed size	Recommended ball charge (pieces)					
			Ø 5 mm	Ø 7 mm	Ø 10 mm	Ø 15 mm	Ø 20 mm	Ø 30 mm
12 ml	up to ≤5 ml	<1 mm	50	15	5	-	-	-
25 ml	up to ≤10 ml	<1 mm	95 – 100	25 – 30	10	-	-	-
50 ml	5 – 20 ml	<3 mm	200	50 – 70	20	7	3 – 4	-
80 ml	10 – 35 ml	<4 mm	250 – 330	70 – 120	30 – 40	12	5	-
125 ml	15 – 50 ml	<4 mm	500	110 – 180	50 – 60	18	7	-
250 ml	25 – 120 ml	<6 mm	1100 – 1200	220 – 350	100 – 120	35 – 45	15	5
500 ml	75 – 220 ml	<10 mm	2000	440 – 700	200 – 230	70	25	8

The table shows the recommended charges (in pieces) of differently sized grinding balls in relation to the grinding jar volume, sample amount and maximum feed size.

PLANETARY BALL MILL PM 200

TYPICAL SAMPLE MATERIALS


RETSCH planetary ball mills are perfectly suitable for size reduction of, for example, alloys, bentonite, bones, carbon fibres, catalysts, cellulose, cement clinker, ceramics, charcoal, chemical products, clay minerals, coal, coke, compost, concrete, electronic scrap, fibres, glass, gypsum, hair, hydroxyapatite, iron ore, kaolin, limestone, metal oxides, minerals, ores, paints and lacquers, paper, pigments, plant materials, polymers, quartz, seeds, semi-precious stones, sewage sludge, slag, soils, tissue, tobacco, waste samples, wood, etc.

Medium-hard: biomass



*35 g sample
125 ml stainless steel grinding jars
7 x 20 mm stainless steel grinding balls
15 min at 500 rpm*

Brittle: tin sulfide

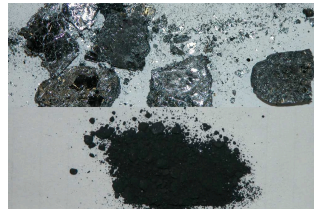

*52 g sample
125 ml agate grinding jars
50 x 10 mm agate grinding balls
60 min at 550 rpm
10 min intervals and direction reversal*

Fibrous: straw



*7 g sample
125 ml zirconium oxide grinding jars
50 x 10 mm zirconium oxide grinding balls
40 min at 300 rpm
10 min intervals and direction reversal*

Hard: natural rocks




*60 g sample
125 ml stainless steel grinding jars
7 x 20 mm stainless steel grinding balls
10 min at 420 rpm*

PLANETARY BALL MILL PM 200

FUNCTIONAL PRINCIPLE

The grinding jars are arranged eccentrically on the sun wheel of the planetary ball mill. The direction of movement of the sun wheel is opposite to that of the grinding jars in the ratio 1:-2. The grinding balls in the grinding jars are subjected to superimposed rotational movements, the so-called Coriolis forces.

The difference in speeds between the balls and grinding jars produces an interaction between frictional and impact forces, which releases high dynamic energies. The interplay between these forces produces the high and very effective degree of size reduction of the planetary ball mill.

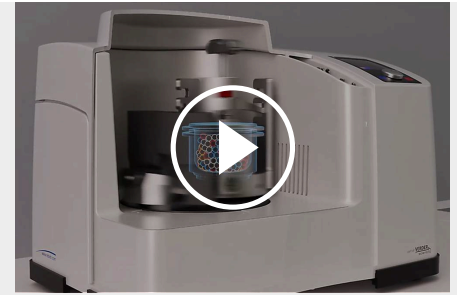

[Click to view video](#)

PLANETARY BALL MILL PM 200

TECHNICAL DATA

Applications	pulverizing, mixing, homogenizing, colloidal milling, mechanical alloying, mechanochemistry, nano grinding
Field of application	agriculture, biology, chemistry, construction materials, engineering / electronics, environment / recycling, geology / metallurgy, glass / ceramics, medicine / pharmaceuticals
Feed material	soft, hard, brittle, fibrous - dry or wet
Size reduction principle	impact, friction
Material feed size*	< 4 mm
Final fineness*	< 1 µm, for colloidal grinding < 0.1 µm
Batch size / feed quantity*	max. 2 x 50 ml
No. of grinding stations	2
Speed ratio	1 : -2
Sun wheel speed	100 - 650 min ⁻¹
Effective sun wheel diameter	157 mm
G-force	37.1 g
Type of grinding jars	EasyFit, optional aeration covers, safety closure devices
Material of grinding tools	hardened steel, stainless steel, tungsten carbide, agate, sintered aluminum oxide, silicon nitride, zirconium oxide
Grinding jar sizes	12 ml / 25 ml / 50 ml / 80 ml / 125 ml
Stackable grinding jars	12 ml / 25 ml
Setting of grinding time	digital, 00:00:01 to 99:59:59
Interval operation	yes, with direction reversal
Interval time	00:00:01 to 99:59:59
Pause time	00:00:01 to 99:59:59
Storable SOPs	10
Measurement of input energy possible	yes
Interface	RS 232 / RS 485
Drive	3-phase asynchronous motor with frequency converter
Drive power	750 W
Electrical supply data	different voltages
Power connection	1-phase
Protection code	IP 30

Power consumption	~ 1250 W (VA)
W x H x D closed	640 x 480 (780) x 420 mm
Net weight	~ 76 kg
Standards	CE
Patent / Utility patent	SafetySlider (DE 202008008473)


*depending on feed material and instrument configuration/settings

www.retsch.com/pm200

ORDER DATA

PLANETARY BALL MILL PM 200

(please order grinding jars and balls separately)

20.640.0001  PM 200 with 2 grinding stations, speed ratio 1 : -2

other electrical versions available for the same price

ACCESSORIES PLANETARY BALL MILLS

22.661.0003  Clamping unit for PM 200

03.025.0178 Adapter for stacking grinding jars 50 ml - 80 ml

02.728.0048  Counter aid for sun wheel PM 100, PM 200 and PM 400

03.486.0062 Opening aid for clamping unit of planetary ball mills

99.200.0008  IQ/OQ Documentation for PM 200

GRINDING JARS EASYFIT

(grinding jars EasyFit are suitable for all planetary ball mills)

HARDENED STAINLESS STEEL

01.462.0239  12 ml

01.462.0240  25 ml

01.462.0516 50 ml

01.462.0517 80 ml

01.462.0518 125 ml

TUNGSTEN CARBIDE

01.462.0494	50 ml
01.462.0495	80 ml
01.462.0527	125 ml

AGATE

01.462.0509	50 ml
01.462.0511	80 ml
01.462.0515	125 ml

SINTERED ALUMINUM OXIDE

01.462.0507	50 ml
01.462.0512	125 ml

ZIRCONIUM OXIDE

01.462.0508	50 ml
01.462.0510	80 ml
01.462.0513	125 ml

ACCESSORIES FOR GRINDING JARS EASYFIT FOR WET GRINDING, GRINDING WITH INERT ATMOSPHERE AND MECHANICAL ALLOYING (MA)

AERATION LIDS (INCL. INLAY)

22.107.0613	for grinding jars EasyFit 50 ml - 125 ml, hardened stainless steel
22.107.0616	for grinding jars EasyFit 50 ml - 125 ml, tungsten carbide
22.107.0617	for grinding jars EasyFit 50 ml - 125 ml, agate
22.107.0615	for grinding jars EasyFit 50 ml - 125 ml, zirconium oxide
22.864.0001	Spare valve set for aeration lids M8x1



INLAY FOR AERATION LID

03.474.0225	for grinding jars EasyFit 50 ml - 125 ml, hardened stainless steel
-------------	--

03.474.0207	for grinding jars EasyFit 50 ml - 125 ml, tungsten carbide
03.474.0208	for grinding jars EasyFit 50 ml - 125 ml, agate
03.474.0206	for grinding jars EasyFit 50 ml - 125 ml, zirconium oxide

AERATION LIDS FOR GRINDING JARS EASYFIT

INCL. O-RINGS AND SINTERED FILTER (PLEASE ORDER LID INSERT AND GRINDING JAR SEPARATELY)

22.107.0636	Aeration lid for grinding jar EasyFit 50 ml - 125 ml
22.107.0637	Aeration lid for grinding jar EasyFit 250 ml - 500 ml

INSERT FOR GRINDING JAR EASYFIT

03.474.0261	Aeration lid insert for grinding jar EasyFit 50, 80 oder 125 ml, stainless steel
03.474.0262	Aeration lid insert for grinding jar EasyFit 50, 80 oder 125 ml, zirconium oxide
03.474.0263	Aeration lid insert for grinding jar EasyFit 50, 80 oder 125 ml, tungsten carbide
03.474.0268	Aeration lid insert for grinding jar EasyFit 50, 80 oder 125 ml, agate
22.186.0007	Sintered filter with O-ring, set of 10 pieces
22.864.0001	Valve set M8x1 for GrindControl and aeration lids



SAFETY CLOSURE DEVICES

22.867.0011	for grinding jars EasyFit 50 ml - 125 ml
02.486.0055	Opening aid for safety closure device

GASKETS FOR GRINDING JARS EASYFIT

O-RINGS

05.114.0086	O-ring for 12 ml grinding jar EasyFit
05.114.0085	O-ring for 25 ml grinding jar EasyFit
05.114.0056	O-ring for 50 ml - 125 ml grinding jars EasyFit
03.111.0438	Flat gasket for 50 ml, 80 ml or 125 ml



GRINDING BALLS

HARDENED STEEL

05.368.0029  5 mm Ø

05.368.0030  7 mm Ø


05.368.0059  10 mm Ø


05.368.0032  12 mm Ø

05.368.0108  15 mm Ø


05.368.0033  20 mm Ø

STAINLESS STEEL

22.455.0010  2 mm Ø, 500 g (approx. 110 ml)

22.455.0011  3 mm Ø, 500 g (approx. 120 ml)

22.455.0002  3 mm Ø, 200 pieces (approx. 6 ml)

22.455.0001  4 mm Ø, 200 pieces (approx. 14 ml)

22.455.0003  5 mm Ø, 200 pieces (approx. 25 ml)

05.368.0034  5 mm Ø

05.368.0035  7 mm Ø

05.368.0063 10 mm Ø



05.368.0037 12 mm Ø



05.368.0109 15 mm Ø



05.368.0062 20 mm Ø



TUNGSTEN CARBIDE

22.455.0006 3 mm Ø, 200 pieces (approx. 6 ml)



22.455.0005 4 mm Ø, 200 pieces (approx. 14 ml)



22.455.0004 5 mm Ø, 200 pieces (approx. 25 ml)



05.368.0038 5 mm Ø



05.368.0039 7 mm Ø



05.368.0071 10 mm Ø



05.368.0041 12 mm Ø



05.368.0110 15 mm Ø



05.368.0070 20 mm Ø




AGATE

05.368.0024 5 mm Ø



05.368.0025		7 mm Ø
05.368.0067		10 mm Ø
05.368.0027		12 mm Ø
05.368.0111		15 mm Ø
05.368.0028		20 mm Ø

SINTERED ALUMINUM OXIDE

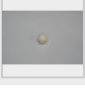05.368.0021		10 mm Ø
05.368.0112		15 mm Ø
05.368.0054		20 mm Ø

ZIRCONIUM OXIDE

32.368.0005		0.1 mm Ø, 0.5 kg (approx. 135 ml)
32.368.0003		0.5 mm Ø, 0.5 kg (approx. 135 ml)
32.368.0004		1 mm Ø, 0.5 kg (approx. 135 ml)
05.368.0089		2 mm Ø, 0.5 kg (approx. 135 ml)
05.368.0090		3 mm Ø, 0.5 kg (approx. 140 ml)
22.455.0007		3 mm Ø, 200 pieces (approx. 6 ml)

22.455.0009



5 mm Ø, 200 pieces (approx. 25 ml)

05.368.0146

7 mm Ø

05.368.0094



10 mm Ø

05.368.0096



12 mm Ø

05.368.0113



15 mm Ø