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The universal standard model

Cutting mills are suitable for the grinding of soft, medium-hard, tough, elastic, fibrous, and heterogeneous mixes of products. Within the group of RETSCH cutting mills, the SM 200 is the universal standard model which covers a vast range of applications.

It is equipped with a strong 2.2 kW motor and operates at a rotor speed of 1.500 rpm. In combination with the optional cyclone-suction-combination, the SM 200 is also suitable for grinding light or heat sensitive sample materials. A heavy-metal-free version is available to avoid cross-contamination.

[Click to view video](#)

FOR UNIVERSAL APPLICATIONS

- | Feed size <60 x 80 mm
- | Defined final fineness due to bottom sieves with aperture sizes from 0.25 - 20 mm
- | Optimum cutting effects thanks to double acting cutting bars
- | Quick and easy cleaning due to smooth surfaces, fold back hopper, push-fit rotor and sieves
- | Operation with gravity outlet or optionally available cyclone unit possible
- | Wide range of accessories including various hoppers, collection systems, rotors and sieves
- | Highest safety standards due to engine brake, central locking device, electronic safety check and safety base frame

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
SUPERIORITY IN DETAIL

DOUBLE ACTING CUTTING BARS

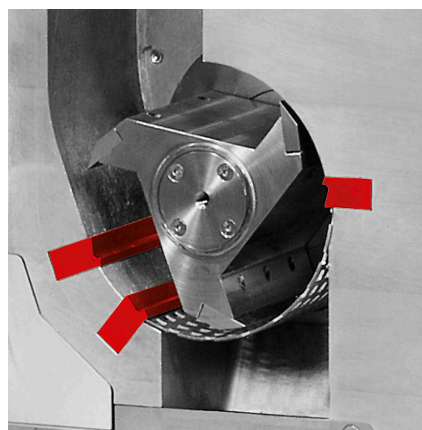Ideal for light & heat-sensitive materials

Instead of the gravity outlet, an optionally available cyclone separator can be connected to the SM 200. The use of the cyclone-suction-combination is advantageous, for example, for light and heat-sensitive samples.

- | Increased throughput
- | Improved sample discharge from the grinding chamber
- | Efficient cooling of sample and cutting tools
- | Collecting vessels with 0.25 - 0.5 - 1 - 2 - 5 and 30 liters are available for the cyclone unit

The SM 200 is equipped with double acting cutting bars which substantially increase the number and effectiveness of the cutting events. The parallel section rotor, for example, generates 18 cutting events per rotation.

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EASY OPERATION & CLEANING

Operation of the SM 200 is exceptionally simple and safe. A safety switch prevents the mill from being switched on with the door open. An electronic safety check ensures that the door cannot be opened when the motor is running – the motor brake makes the rotor come to an immediate standstill after the mill has been switched off.

The hopper can be folded back to get full access to the grinding chamber which, in combination with the smooth surfaces, greatly facilitates cleaning the mill. The push-fit rotor can be removed without tools.

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ACCESSORIES FOR SAFE AND EFFECTIVE GRINDING PROCESSES

The SM 200 can be equipped with different rotors, hoppers and bottom sieves. The mill and accessories, including cutting bars and cutting plates, are also available in versions for heavy-metal-free size reduction.


ROTOR TYPES

A parallel section rotor and a 6-disc rotor are available in stainless and heavy-metal-free steel for optimum adaptation to sample properties.

HOPPER TYPES

The SM 200 can be equipped either with the universal hopper, which is suitable for most samples, or the long-stock hopper which is specially designed for long samples.

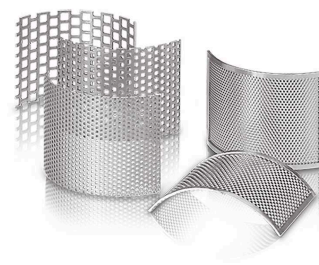
BOTTOM SIEVES

The bottom sieves with trapezoid or square holes are made of stainless steel or heavy-metal-free steel.

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
RETSCH cutting mills are suitable for a vast range of applications. Typical materials include PET preforms, aluminium slag, animal feed, bones, cables, cardboard, electronic components, feed pellets, foils, food, leather, lignite, material mixtures, non-ferrous metals, paper, pharmaceutical products, plant materials, plastic toys, plastics, polymers, refuse derived fuels, resins, rubber, spices, straw, textiles, waste, wood, etc.

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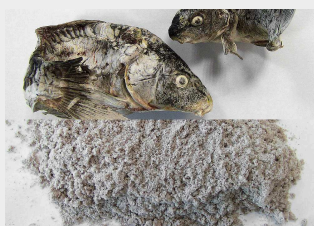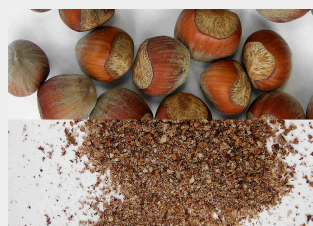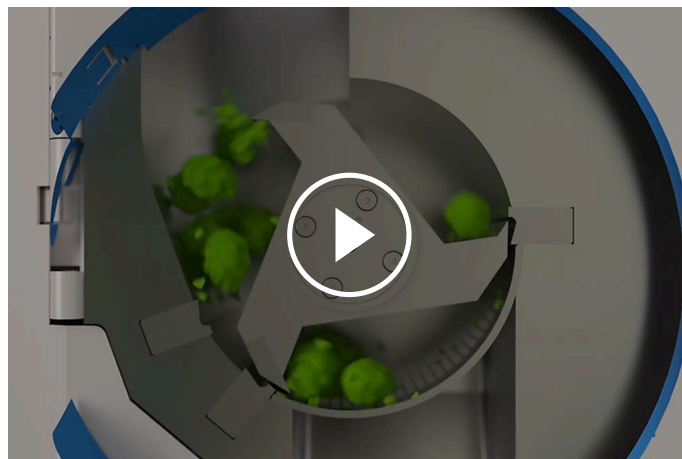
CUTTING MILL SM 200

Cutting Mill SM 200 complete with collecting receptacle, 5 litres and base frame with wheels
(please order rotor, hopper and bottom sieve separately)

20.728.0001		SM 200	3/N~ 400 V, 50 Hz
20.728.0002		SM 200	3/N~ 400 V, 50 Hz, for grinding without heavy-metal contamination

other electrical versions available for the same price

ROTORS SM 100, SM 200, SM 300

22.608.0022		6-disc rotor, stainless steel, with reversible cutting tips of tungsten carbide
22.608.0028		6-disc rotor, steel 1.0580, with reversible cutting tips of tungsten carbide, for grinding without heavy-metal contamination
22.608.0021		Parallel section rotor, stainless steel
22.608.0018		Parallel section rotor, steel 1.0580, for grinding without heavy-metal contamination

HOPPERS

22.785.0001		Universal hopper with plastic plunger for SM 200 / SM 300
22.408.0003		Long stock hopper with wooden plunger for SM 200 / SM 300

ADDITIONAL PLUNGERS FOR UNIVERSAL HOPPER

22.725.0010		Plastic plunger
22.725.0011		Wooden plunger
22.725.0012		Aluminum plunger

BOTTOM SIEVES SM 100, SM 200, SM 300

BOTTOM SIEVES OF STAINLESS STEEL

03.647.0312		Trapezoid holes	0.25 mm
03.647.0313		Trapezoid holes	0.50 mm
03.647.0314		Trapezoid holes	0.75 mm
03.647.0315		Trapezoid holes	1.00 mm
03.647.0316		Trapezoid holes	1.50 mm
03.647.0317		Trapezoid holes	2.00 mm
03.647.0318		Square holes	2.00 mm
03.647.0319		Square holes	3.00 mm
03.647.0320		Square holes	4.00 mm
03.647.0321		Square holes	6.00 mm
03.647.0322		Square holes	8.00 mm
03.647.0024		Square holes	10.00 mm
03.647.0062		Square holes	20.00 mm

BOTTOM SIEVES OF STEEL 1.0353 FOR GRINDING WITHOUT HEAVY-METAL CONTAMINATION

03.647.0325		Trapezoid holes	0.25 mm
-------------	---	-----------------	---------

03.647.0326		Trapezoid holes	0.50 mm
03.647.0327		Trapezoid holes	0.75 mm
03.647.0328		Trapezoid holes	1.00 mm
03.647.0329		Trapezoid holes	1.50 mm
03.647.0331		Square holes	2.00 mm
03.647.0332		Square holes	3.00 mm
03.647.0333		Square holes	4.00 mm
03.647.0334		Square holes	6.00 mm
03.647.0335		Square holes	8.00 mm
03.647.0359		Square holes	10.00 mm

ACCESSORIES SM 100, SM 200, SM 300

CYCLONE-SUCTION-COMBINATION FOR SM 200 / SM 300

22.020.0004	 	Cyclone unit incl. 1 sample bottle 500 ml (to be used with industrial vacuum cleaner)
22.935.0039		Cyclone unit incl. 30 litres collection receptacle (to be used with industrial vacuum cleaner). Requires 02.183.0156 cyclone outlet.
22.748.0005		Industrial vacuum cleaner for mills, 230 V, 50/60 Hz (other electrical versions and spare filter bags available upon request)
22.523.0002		Sample bottles, 500 ml, 10 pieces
22.001.0001		Holder for sample bottles 1 / 2 / 5 litre
05.239.0011		Sample bottle, 1 litre, 1 piece, to be used with holder 22.001.0001
05.239.0010		Sample bottle, 2 litre, 1 piece, to be used with holder 22.001.0001

05.239.0009  Sample bottle, 5 litre, 1 piece, to be used with holder 22.001.0001

22.003.0014  Collecting receptacle, 5 litre, stainless steel, incl. adapter lid

COLLECTING RECEPTACLES

22.003.0008  Collecting unit 250 / 500 ml, incl. 2 sample bottles

22.523.0001  Sample bottles, 250 ml, 10 pieces

22.003.0007  Collecting receptacle, 5 litres, stainless steel

05.045.0528 Collecting receptacle, 30 litres

22.187.0006   Ring-type filter for collecting receptacle, 5 litres

22.748.0001  Dust filter clamping rings for ring-type filter, incl. 5 dust filters

22.524.0002   Dust filters for ring-type filter, 25 pcs.

22.187.0007  Filter hose for collecting receptacle, 5 litres, incl. holder

ACCESSORIES FOR ROTORS

22.908.0001  Reversible cutting tip incl. screws, tungsten carbide, for 6-disc rotor (10 pcs.)

22.151.0004  Cutting plates, stainless steel, for parallel section rotor, 1 set (3 pieces)

22.151.0006   Cutting plates, steel 1.1740, for parallel section rotor, 1 set (3 pieces) for grinding without heavy-metal contamination
for grinding without heavy-metal contamination

STATIONARY CUTTING BARS

22.152.0003  Double acting cutting bars, stainless steel, 1 set (3 pieces) for SM 200 / SM 300

22.152.0005   Double acting cutting bars, steel 1.1740,
1 set (3 pieces)
for SM 200 / SM 300 for grinding without heavy-metal contamination

22.608.0035 Adjustment aid for cutting bars

ADDITIONAL ITEMS SM 100, SM 200, SM 300

22.225.0002  Replacement extraction tool for rotor

22.111.0003 Gasket for door SM 100 / SM 200 / SM 300, 1 piece

99.200.0013  IQ/OQ Documentation for SM 200

99.200.0014  IQ/OQ Documentation for SM 300